

# Palliative Care Tool User Guide V1.2

The Palliative Care Toolbar supports clinicians in earlier identification of patients with a life-limiting illness who could benefit from a palliative care approach



## **Table of Contents**

Intended Tool Use	
1. Version History	2
2. Toolkit Contents	3
Installation	5
Installation process	5
Importing Custom Forms and Encounter Assistants	6
2. Installing searches	8
3. Setting a Reminder to show the Palliative Toolbar	10
4. Importing Handouts into PSS	12
User Guide	13
The "No" pathway for the Palliative Tool	13
The "Yes" pathway for the Palliative Tool	16
The "Not Sure" pathway for the Palliative Tool	18
How to use the Palliative – Main Form	22
What is in the View: panel?	22
Supporting tools in the Palliative – Main Form	26
What is the Palliative SOAP form?	26
Other helpful supporting tools in the Palliative – Main Form	27
Updating/Changing resources	28

## **Intended Tool Use**

The Palliative Care Toolbar supports clinicians in earlier identification of patients with a life-limiting illness who could benefit from a palliative care approach. This tool can provide decision support in assessing the palliative needs of the patient and is then provided with the appropriate resources, relative to the patient's current Palliative Performance Scale (PPS) phase. The tool is also equipped with Advance Care Planning tools that allow clinicians to customize a reassessment date to prompt the Palliative Care tool to be used again.

**Important note:** this tool was originally created by the eHealth Centre of Excellence, which rebranded to Amplify Care in 2025.



# **1. Version History**

Version	Release Date	Updates	
1.2	September 2025	<ul> <li>Search criteria now looks at history of past health in the CPP</li> <li>Resources have been updated/refreshed</li> <li>Hidden resource library to customize with pre-configured resources now included</li> <li>Ability to update/customize resource list remains</li> <li>Resources organized by patient-facing, clinician facing and regionally-specific resources</li> </ul>	
1.0	2018	Originally created from a partnership between representatives from the Waterloo Wellington Integrated Hospice Palliative Care Regional Program and Amplify Care (formerly known as The eHealth Centre of Excellence)	



#### 2. Toolkit Contents

The palliative care toolkit includes:

#### **Custom forms**

- Advance Care Planning Template.cfm
- Palliative Criteria Match.cfm
- Palliative BPI.cfm
- Palliative General Indicators.cfm
- Palliative Main Form.cfm
- Palliative CHF Indicators.cfm
- Palliative CKD Indicators.cfm
- Palliative MRC Indicators.cfm
- Palliative Office Visits Reference.cfm
- Palliative PHQ9.cfm
- Palliative Reassess.cfm
- Palliative SOAP.cfm
- Palliative Special Premium Bonus.cfm
- Palliative Specific Indicators Reference.cfm
- Palliative Specific Indicators.cfm
- Palliative Toolbar.cfm
- Palliative Performance Scale.cfm
- Substitute Decision Maker SDM.cfm
- Symptom Response Kit Order Form.cfm
- Palliative ESAS-R.cfm
- Palliative PPS Instructions.cfm
- Palliative General indicators reference.cfm
- Amplify-FeedbackForm.cfm
- Amplify-InitialTermsForm.cfm
- Amplify-MainForm-OptedIn.cfm
- Amplify-MainForm-OptedOut.cfm
- Amplify-NewVersionForm.cfm

#### <u>Handouts</u>

- Palliativecarebillingguide.pdf
- DNRC Information Bulletin

#### Resource Library handouts

#### Reminders and Searches

- Palliative Toolbar Reminder.srx
- Palliative All Flagged Patients.srx
- Palliative Identified Patients.srx
- Palliative-Criteria-Age.srx
- Palliative-Criteria-ALS.srx
- Palliative-Criteria-ALS-HPH.srx
- Palliative-Criteria-Cancer.srx
- Palliative-Criteria-Cancer-HPH.srx
- Palliative-Criteria-CHF.srx
- Palliative-Criteria-CHF-HPH.srx
- Palliative-Criteria-CKD.srx
- Palliative-Criteria-CKD-HPH.srx
- Palliative-Criteria-COPD.srx
- Palliative-Criteria-COPD-HPH.srx
- Palliative-Criteria-Dementia.srx
- Palliative-Criteria-Dementia-HPH.srx
- Palliative-Criteria-Dyspagia.srx
- Palliative-Criteria-Dyspagia-HPH.srx
- Palliative-Criteria-Frailty.srx
- Palliative-Criteria-Frailty-HPH.srx
- Palliative-Criteria-LiverDisease.srx
- Palliative-Criteria-LiverDisease-HPH.srx
- Palliative-Criteria-MS.srx
- Palliative-Criteria-MS-HPH.srx
- Palliative-Criteria-NYHA3or4.srx
- Palliative-Criteria-Palliative.srx
- Palliative-Criteria-Palliative-HPH.srx
- Palliative-Criteria-Parkinsons.srx
- Palliative-Criteria-Parkinsons-HPH.sr
- 7 Things to know about palliative care
- 10 myths about palliative care



Install these handouts if you'd like to use the optional resource library to include them in your tool. You can preview them by opening the files found in the zipped toolkit folder.


Users can also add their own resources (websites, custom forms and handouts) using the built in instructions shown when opening the main form in editor mode. \*Note you may need a custom form license\*

- · Hospice Room pamphlet
- maid-individuals-and-families
- MEDICAL ASSISTANCE IN DYING (MAID)
- ontario caregiver starter kit
- · serious illness conversation guide
- SPIRITUAL ASSESSMENT TOOL- FICA

# Installation


## Installation process

1. Select the Palliative Toolkit zip file attachment and save it to your desktop or another file folder location of your preference.

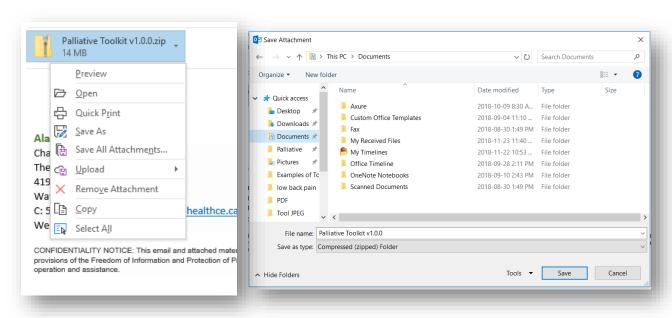


Extract and save the Palliative Toolkit zip file to your desktop or other file folder location
of preference. Right-click on the zip file and select Extract All... to save to the file
location.



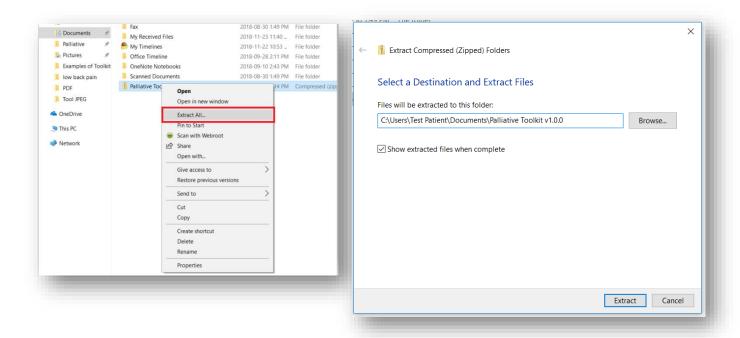


## 1.Importing Custom Forms and Encounter Assistants

1. Login to Practice Solutions and select the **Records** Tab

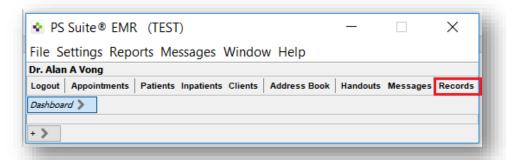





2. In the Records window, select the Settings tab -> Edit Custom Forms

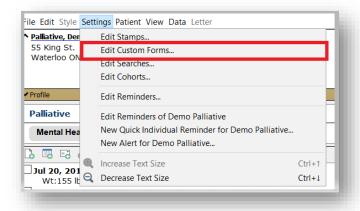

3. Select the File tab -> Import Form(s)...

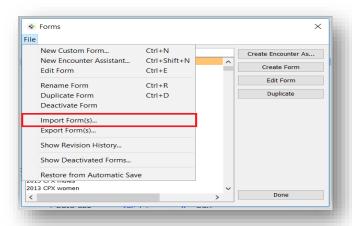


3. Find the file location where you have saved the unzipped Palliative Toolkit files, select all by "Ctrl + A" or "Apple button + A" so it automatically Select all. Click **Choose**.

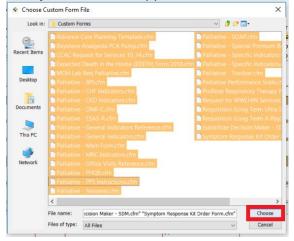




- 4. To insert one of Palliative files into a patient's chart/record (i.e. the Palliative Main Form):
  - a. Find the patient record in PSS

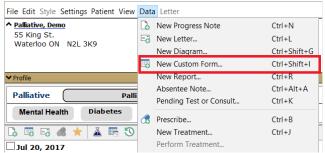

- b. Select the Data tab -> New Custom Form...
- Search for the custom form you would like to insert (i.e start typing keyword
   Palliative Main Form in the text field) and select it from the custom form list ->
   Choose This Form.

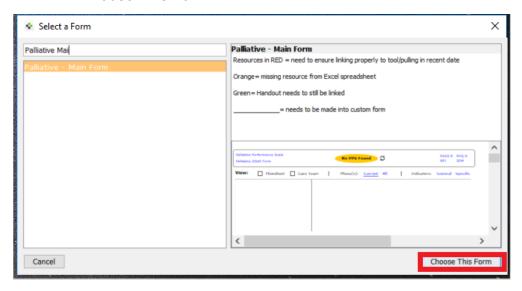


## 2. Installing searches


**Notes - 1.** It is recommended to install the search files ahead of setting up the toolbar reminder as the toolbar reminder references many of the search files and installing out of order may trigger an error. **2.** If you are updating from an older version — It is recommended to delete all existing versions of the 'palliative - ###' searches to reduce duplicates being added to your EMR system.

1. From the Records window, choose Settings > Edit Searches.



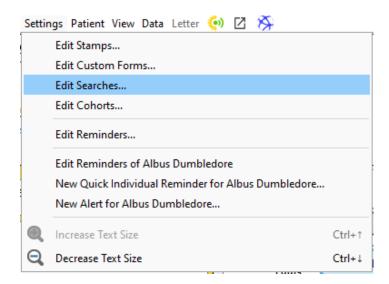


2. in the Edit Searches window, from the Edit menu, choose Import Searches.

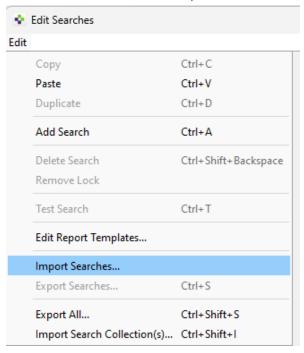


3. Find the file location where you have saved the unzipped Palliative Toolkit files, select all by "Ctrl + A" or "Apple button + A" so it automatically Select all. Click **Choose**.



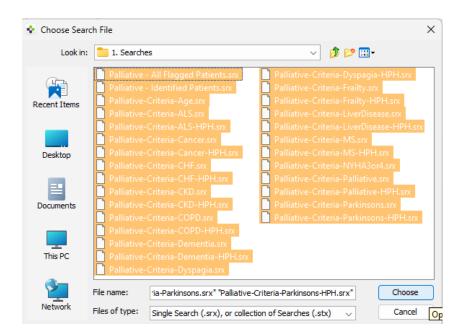


#### 3. Setting a Reminder to show the Palliative Toolbar

1. Log in to Practice Solutions and select the Records tab.

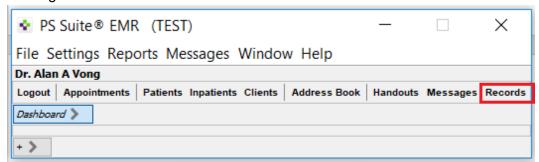


2. Select the Settings tab -> Edit Reminders...

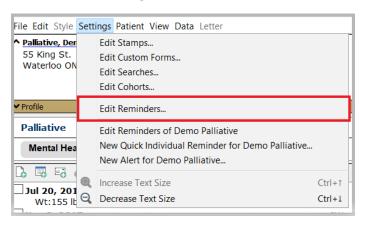





3. Select the File tab -> Import Form(s)...

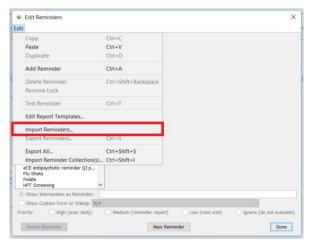

4. Find the file location where you have saved the unzipped Palliative Toolkit files, select all by clicking on one file and holding Shift and selecting the others in the list. Click **Choose**.

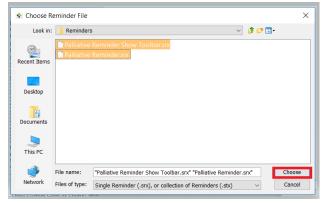


5. Search for Palliative Reminder Show Toolbar, then select Custom Form or Stamp and find the **Palliative – Toolbar** custom form in the dropdown. Click done.

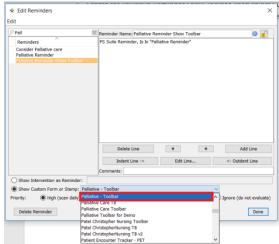




## 4. Importing Handouts into PSS


1. Log in to Practice Solutions and go to → Handouts

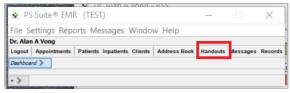

2. Select Edit → Import Handout

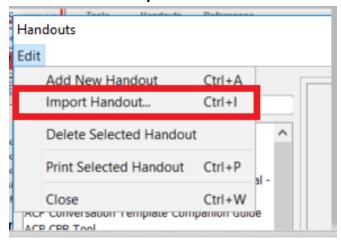

Locate the Palliative Toolkit v1.0.0 folder and find the Handouts. Have the Window with the handouts open and drag each handout into the Handout section (can only drag one at a time). Click **Done**.

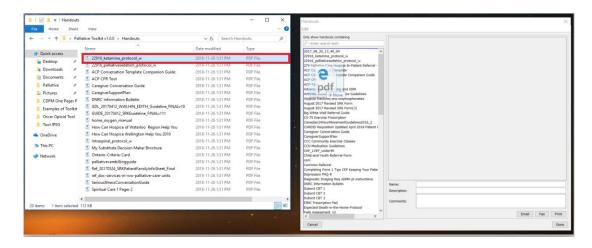




# **User Guide**

# The "No" pathway for the Palliative Tool

1. Once everything has been installed, the Palliative toolbar should appear in the patient's chart. Click on the button called "Palliative?"




2. The Surprise Question will prompt as a pop-up and the options: **No**, **Not Sure**, **Yes and Cancel** will appear. Click **No** if you are not surprised about your patient.

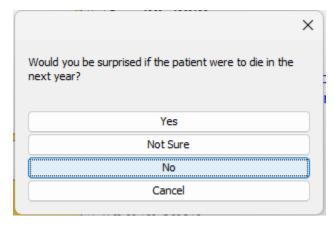




 A pop-up will appear when the Palliative Performance Scale has not been completed. The completion of Palliative Performance Scale will determine the appropriate patient's phase and resources.

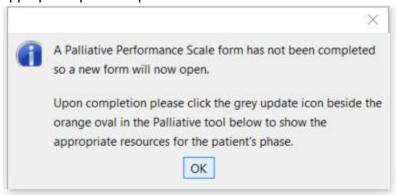


4. Fill in the Palliative Performance Scale (PSS) and click Add to Notes.

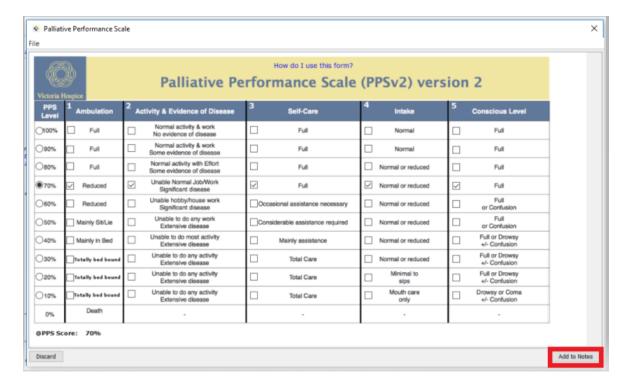



5. Look at the **Palliative – Main Form** and click the refresh button to update the custom form. The **Palliative – Main Form** will pull the latest **PSS** score.

#### **Before:**

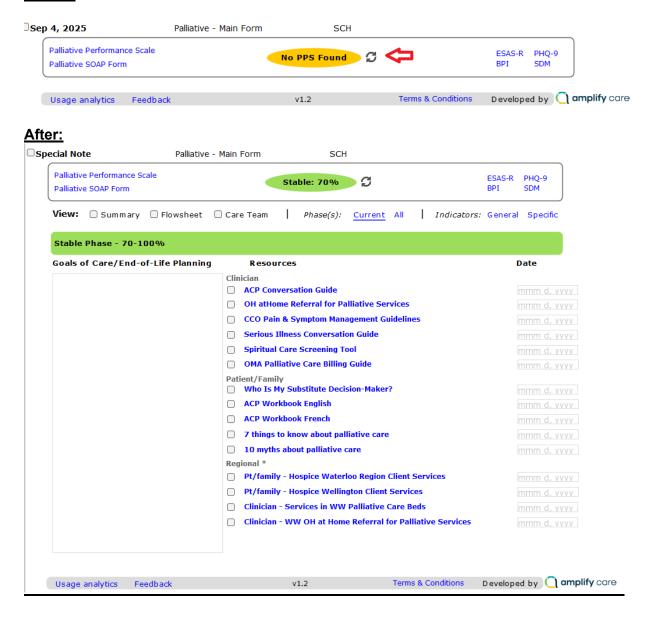




# The "Yes" pathway for the Palliative Tool

1. Once everything has been installed, the Palliative toolbar should appear in the patient's chart. Click on the button called "Palliative?"




2. The Surprise Question will prompt as a pop-up and the options: **No**, **Not Sure**, **Yes and cancel** will appear. Pick and choose **Yes** if you are surprised about your patient.

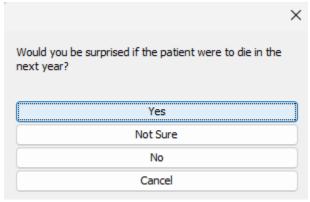





3. A custom form will appear called **Palliative – Reassess.** The **Advance Care Planning** form allows the patient to be reassessed in an appropriate timeframe. Click **Finish**.

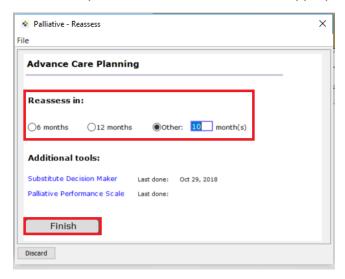

4. In the toolbar section, it will display a reminder button and count down till reassessment is due.

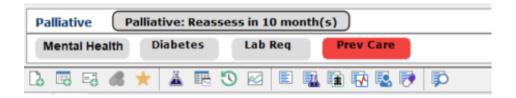

5. The button will turn yellow when it is time for palliative reassessment. Click the button to answer the Surprise Question again.

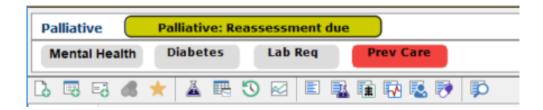




# The "Not Sure" pathway for the Palliative Tool

1. Once everything has been installed, the Palliative toolbar should appear in the patient's chart. Click on the button called "Palliative?"




2. The Surprise Question will prompt as a pop-up and the options: **No**, **Not Sure**, **Yes and cancel** will appear. Click **Not Sure** if you are not sure about your patient.

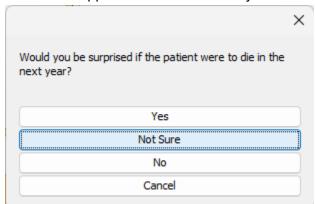





3. A custom form called **Palliative – General Indicators** will display. Click **Yes** if there are general indicators of decline and increasing needs. Click **No** otherwise.

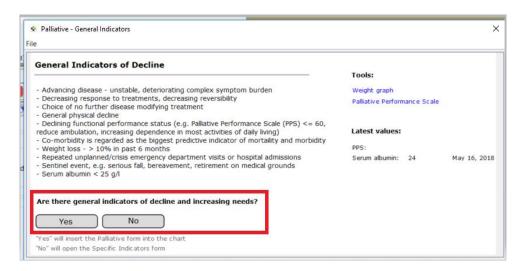

4. If Yes is clicked, a Palliative Performance Scale will need to be filled.

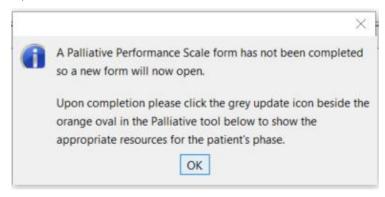

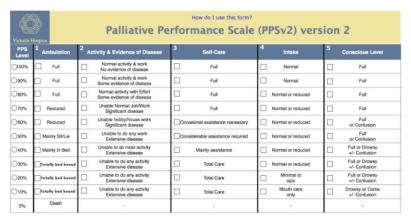

5. If **No** is clicked, a custom form called **Palliative – Specific Indicators** will display. This form provides specific indicators to review about patient.







6. Click **Yes** or **No** once determining the appropriate next step. "**Yes**" will insert the Palliative form into the chart and "**No**" will open the Advance Care Planning form.






7. If **Yes** is clicked, a **Palliative Performance Scale** will need to be filled.

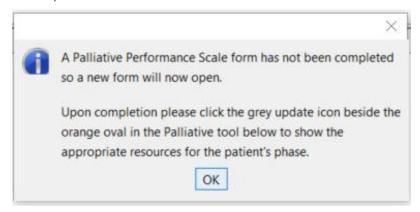

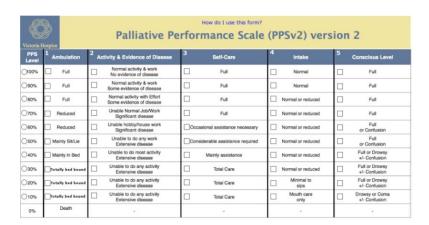

If No is clicked, a custom form will appear called Palliative – Reassess. The Advance
 Care Planning form allows the patient to be reassessed in an appropriate timeframe.
 Click Finish.

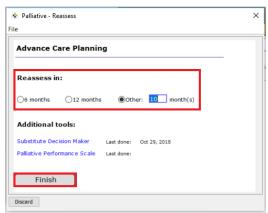



## How to use the Palliative - Main Form


#### What is in the View: panel?

The **View:** panel will contain the core sections of the Palliative tool, sections with checkboxes can be unchecked to minimize space. Blue text is clickable and will direct to one of the following: website, PDF or custom form.

1. Click **Flowsheet** to view the values that are being tracked throughout time by the **Palliative – Main Form**.




Click Care Team to view and document current providers on the patient's care team (the Other fields can be edited when clicked).

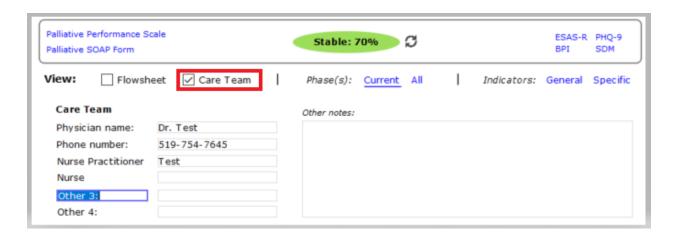




Click <u>Current</u> to view the appropriate phase for the patient (based on the <u>Palliative</u> <u>Performance Scale</u>).

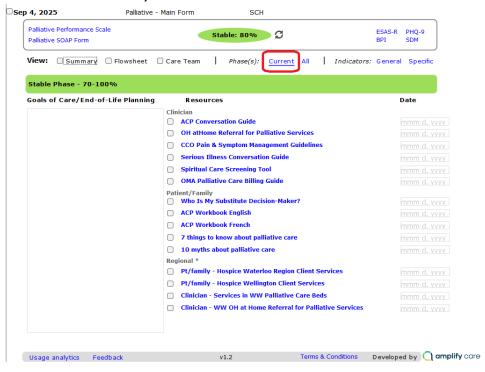

Under **Resources**, click on checkboxes to document the resources used with the patient. It will log the date under **Date**.

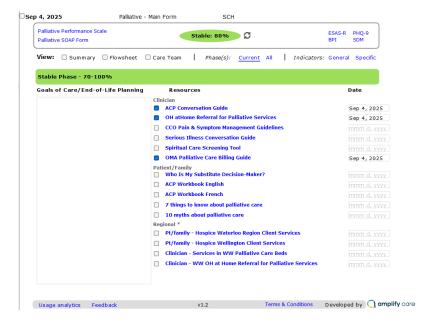




4. Click **All** to view all the palliative care phases of the patient.




5. Click **Copy goals from previous stage** to copy from a previous stage. The same **Resources** can be in different stages and will be documented if selected.

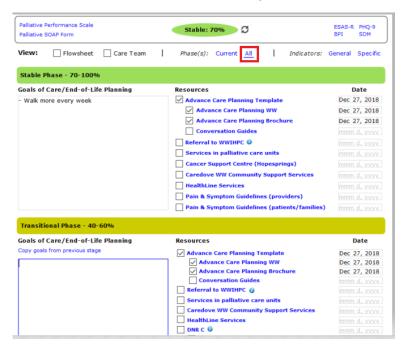




6. Click **General** to view the **Palliative – General Indicators Reference** custom form. It will provide information about the general indicators of decline.

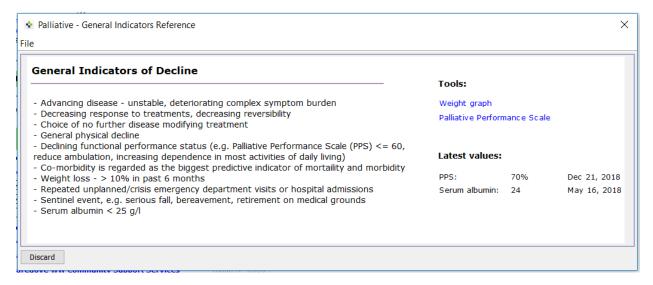

7. Click **Specific** to view the **Palliative – Specific Indicators of Decline** custom form. It will provide information about the specific indicators of decline.

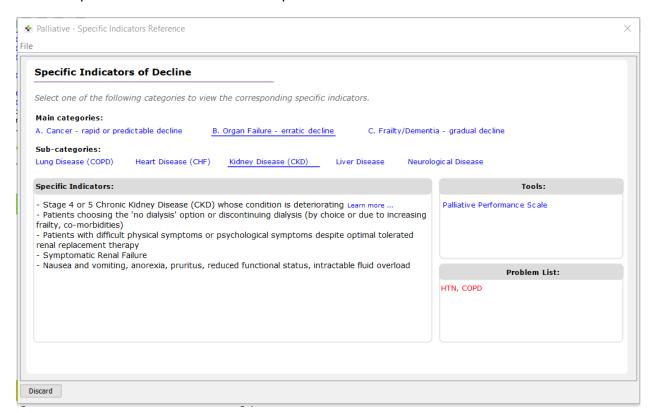




# Supporting tools in the Palliative – Main Form

#### What is the Palliative SOAP form?


1. Click Palliative SOAP Form to view eligible billing codes and document in SOAP notes.



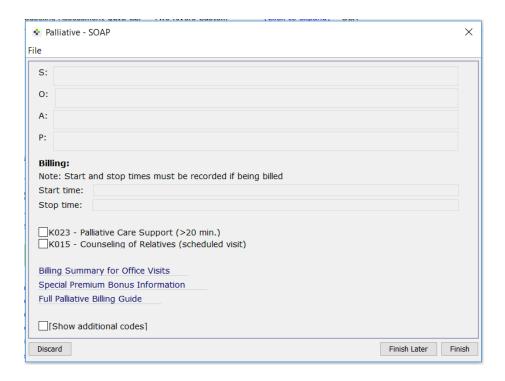



2. Click **Show additional codes** to view all the common forms attached with billing codes.

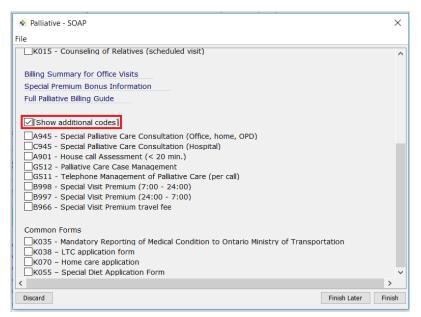


## Other helpful supporting tools in the Palliative - Main Form

1. Click the highlighted box for resources that helps support palliative care.





2. Click **SDM.** This is an example of a resource to use with patients. Fill out the patient's substitute decision maker information. Click **Add to Notes.** 

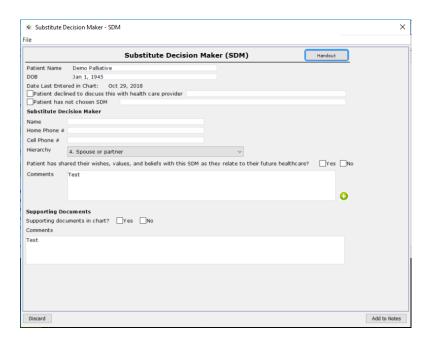

# **Updating/Changing resources**

Resources in the palliative main form are organized by Clinician, Patient/family, and Regional categories. Resources in the Clinician and Patient/Family categories are more generalized and applicable to most regions. Items in the Regional category can be tailored specific to regional care pathways and resources. If you would like to adjust these resources you can:

\*please note you or someone at your site may need to have access and knowledge in custom form editing. All edits must be done in the custom form editing mode of your EMR.


1. Use the built in instructions to add new resources – instructions can be found when opening the form in the custom form editor mode



Add Your Own Resources Guide					
Samples Resources:					
Handout: Web Link: Custom Form:	Conversation Guides Advance Care Planning WW ACP Template	For each of the sample resources, look at the `Click Event` and `Resource` properties as a guide to configure your own resource			
Checkbox(es):  These checkboxes have been setup but still need the `Resource Name` and `Name` to be configured					
< Pre-configured for stable phase					
< Pre-configured for transitional phase					
Pre-configured for end-of-life phase					
Choose a common `Resource Name` for each checkbox. E.g. MyResource     For each checkbox, set the `Name` field in this specific way:      chk + (Resource Name) + (phase)					
e.g. chkMyResourceS					
Phase naming: S = stable, T = transitional, E = End-of-life					
3. The `Event Trigged` event should be `processCheck` for each checkbox					
Date Field(s): These date fields have been setup and only need the Name to be configured					
< Pre-configured for stable phase					
< Pre-configured for transitional phase					
	< Pre-configured for end-of	-life phase			
1. For each date field, set the `Name` field in this specific way:					
date + (Resource Name) + (phase)					
e.g. dateMyResourceS					
Phase naming: S = stable, T = transitional, E = End-of-life					

2. Use our pre-built resource library – items in this section can be swapped into applicable sections of the main form.



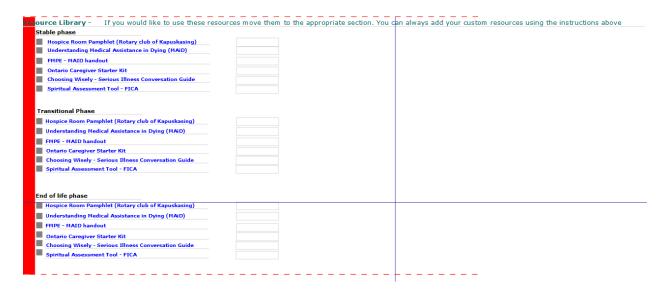

If you have questions or feedback, please feel free to reach out to <a href="mailto:EMRtools@amplifycare.com">EMRtools@amplifycare.com</a>.